



Easy 4 Step Installation

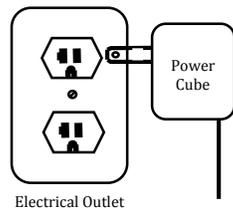


Step 1.

Plug the Power Cube into an outlet.

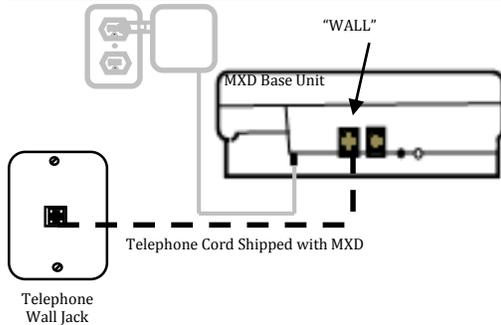


Step 2.

Take the loose end of the gray phone cord coming from the back of the unit and plug it into wall jack. (If your telephone cord is already plugged into the wall jack, remove it first.)

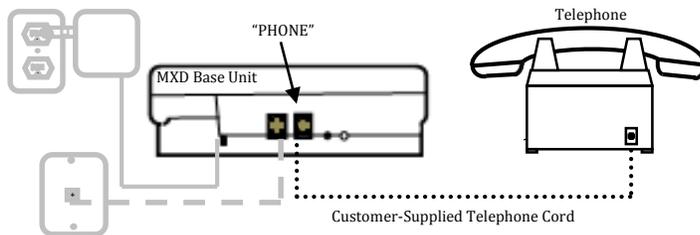


Step 3. (Optional)

Take phone cord that was removed from wall jack (Step 2) and plug into back of the unit marked PHONE.

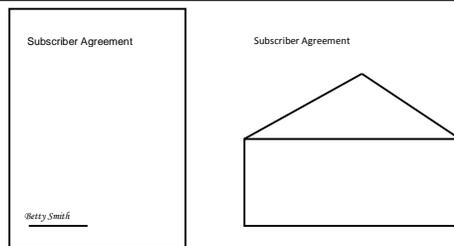

Step 4.

Complete Subscriber Agreement and mail to:

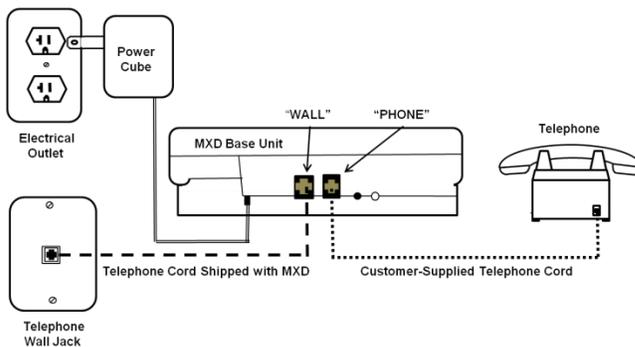

Completed Installation

After installation is complete, press your button on the necklace to do a live test. In about 30 to 45 seconds, a care specialist will be heard over the unit's speaker to welcome you.